

Air Boots Massager

Instruction Manual

ITEM: BH003F



Air Boots Massager operating manual read all INSTRUCTIONS BEFORE USE.

SPECIFICATION

Input power : AC220V / 50Hz

Output power : Dc 7.5V

Electric current : 500mA

Size : 50 x 25 x 3cm

Weight : 1.4Kg

Air pressure higher than 600mmHG

If accident happened because of the human-made, our company couldn't take on any responsibility. This product will be happened the heavy accident because of the wrong usage, therefore please follow this user manual. Use this machine correctly upon the responsibility.


This manual describes the 2 dangerous level, please understand adequately and follow this user manual to use this product.

Attention

1. This product is the foot massager for home use, it COULD NOT for business. And this product could ONLY use on your foot.
2. When you wear this product, one set for one time. Because if your wear this product over 1 set, it will be dangerous.
3. The guardian should be notice especially, the children couldn't use this product while no adult be with them. Also, this product is NOT a toy. And don't let any animal to approach this product.
4. This product could NOT be used over 1 person one time.
5. When you wear the pacemaker, the artificial heart-lung and ECG machine, please DO NOT use this product, it will lead to the mistake.
6. While you use this product, please wear it in our leg. If you wear it in your head or chest, you might suffer serious damage.
7. Do not use it when you would like to sleep or you're sleeping.
8. The strong massager, it is not only ineffective but also your muscle will be ache. So, please adjust the suitable intensity to use. When you feel it is too strong, you could pull the pipe out to reduce the air press intensity.
9. If you would like to use it over 20 minutes, please pull the adaptor from the outlet. And use it after 30 minutes.
10. The person who are taking the treatment or the following person should be follow the doctor instruction to use it:
 - (1) The human who is with spinal curve or spinal abnormality.
 - (2) The human who has bad heart or the heart has malignant tumor.
 - (3) The human who is in pregnancy or in menstrual.
 - (4) The human with skin disease
 - (5) The human who with high body temperature and hypertension
 - (6) The human who is with fracture, dislocation and muscle pain.
 - (7) The human who need quite rest.
 - (8) The human who is with acute irritation and swollen pain.
 - (9) The human who is has heart disease and myocardial infarction.
 - (10) The human who is has perceptive obstacle.
 - (11) The human who feel not good.
11. When you use it, you are not comfortable, vertigo or heartbeat fast. Please stop using immediately and go to see the doctor.

Trouble shooting

- 1.No power, please check if the adaptor is pulled and broken.Also, if using our adaptor.
2. Have the movement but don't feel the air pressure: Open the power, it needs some time to feel air pressure. Please check the condition for minutes. Also, check if the air bag is loaded correctly. If you don't wear the air bag correctly, the air pressure will be focused on one air bag. Or if the air bag is pressed down or stepped or snapped or damaged.

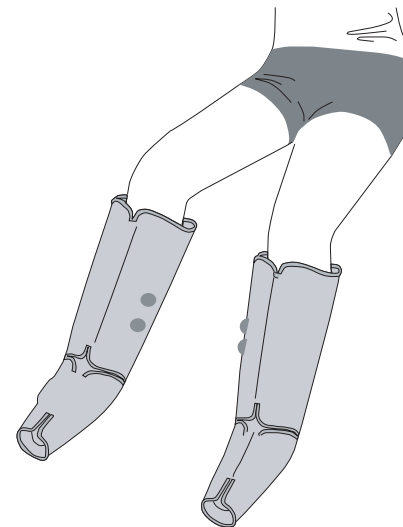

Description on functions of air pressure pants controller

1. The air pressure bag of this controller are controlled separately .
2. Pressing the air pressure switch will start air pressure pants of which there are three modes for adjustment,e.g.,Weak,Mid,Strong,Off.
3. If the switch is set at Weak, the green indicating light will come on. After the air pressure pants is stopped for 15 seconds,it will be inflated automatically for 20 seconds,and so on.After operation for 20 minutes,it will be shut off automatically.
4. If the switch is set at Mid,the yellow indicating light will come on. After the air pressure pants is stopped for 15 seconds,it will be inflated automatically for 30 seconds,and soon.After operationfor 20 minutes,it will be shut off automatically.
5. If the switch is set at Strong,the red indidating light will come on. After the air pressure pants is stopped for 15 seconds,it will be inflated automatically for 37 seconds,and so on.After operationfor 20 minutes,it will be shut off automatically.
6. It is recommended that during initial use, the air pressure pants should be adjusted to Weak, after it is used for a period of time,it can be adjusted to Mid or Strong depending on the condition of the user.
7. The air pressure of this poduct are controlled separately so that the user can start either one depending on his(her) own requirement.
8. It is recommended that this product shall be used for 50 minutes a day and after use,please drink 500c.c. water thus helping body slimming significantly.

Strong Mid,Weak air pressure indicating lights

